

BetterLiving

Wheeled Walker Budget



Soft touch braking mechanism

Durable moulded orthopedic handgrip

Height adjustable handles (795-950mm)

Foam backrest

Steel frame with tube-in-tube construction for increased strength

Vinyl covered foam seat with concealed PVC seat pouch

Removable basket

Compact fold down design

153 mm (6 inch) castors; double locking rear assembly for greater durability and stability

The Budget Wheeled Walker is the economical choice for general indoor use. The sturdy steel frame increases product strength and durability while six inch castors ensure easy operation and manoeuvrability on interior surfaces.

The durable moulded handgrips with soft touch brakes are ideal for people with limited strength or mobility; maximising user comfort by minimising pressure on the hand. Handle height can be adjusted easily to best suit the user, with an extension range of 155 mm to suit users of average weight and height.

BetterLiving

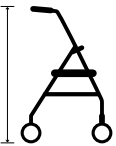

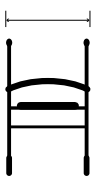
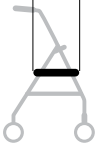
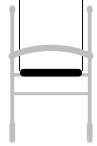
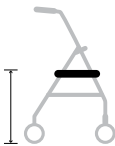


Wheeled Walker Budget



The high quality vinyl covered foam seat has no sharp edges, reducing the potential for skin tear, plus foam backrest for added comfort when seated. Cleverly concealed beneath the seat is a PVC pouch, offering additional security for small personal valuables. The under seat basket is also removable, for the security and transport of goods when leaving the wheeled walker unattended.

Easy closing mechanism and compact fold down design for convenient transport or storage

Colour	Metallic Black
Code	BL7388
Capacity	130 kg
Weight	9 kg
Conforms to	ARTG 176830 ISO 11199-3:2005
Warranty	1 year

HEIGHT	LENGTH	WIDTH	SEAT DEPTH	SEAT WIDTH	SEAT HEIGHT	FOLDED WIDTH	CASTOR SIZE
							
950 mm max 795 mm min	620 mm	610 mm	320 mm	360 mm	545 mm	220 mm	153 mm 6 inch

THERAPEUTIC DEVICES SHOULD ONLY BE USED IN ACCORDANCE WITH MANUFACTURERS INSTRUCTIONS AND UNDER THE CONSENT, SUPERVISION AND MANAGEMENT OF A SUITABLY QUALIFIED HEALTH PROFESSIONAL.