

The Omron Standard **Blood Pressure** Monitor is a compact, fully automatic upper arm blood pressure monitor that uses the oscillometric principle to simply and accurately measure systolic and diastolic blood pressure and pulse rate.


With an easy fitting arm cuff and single button operation, the compact unit offers a simple and convenient way to monitor blood pressure at home. With the touch of a button, the lightweight desktop monitor simultaneously displays systolic, diastolic and pulse rate.

'Enhanced IntelliSense Technology' automatically applies the right amount of pressure for fast, accurate and more comfortable measurements.

A 30 memory storage capacity allows the review of the last 30 readings (including time and date).

Also features an irregular heartbeat detection, to alert if systolic or diastolic pressure is outside the standard range (above systolic 135mmHg/diastolic 85mmHg) plus hypertension indicator and body movement detection.

OMRON

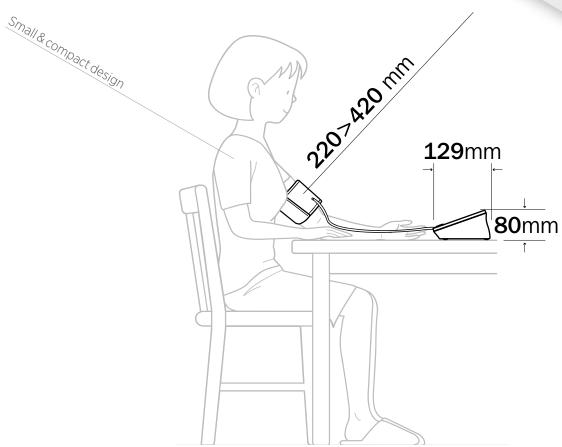

- Large cuff with cuff wrapping guide included.
- Requires 4 x AA batteries (included)
- Measures up to 1,000 times (with new batteries).
- AC power adaptor optional
- Reading time is approximately 50 seconds.

CONSULT YOUR MEDICAL PROFESSIONAL TO EVALUATE THE READINGS. CHECK YOUR

ITEM	Blood Pressure Monitor	AC Power Adaptor
CODE	HEM7121-Standard	BP-Pack
CUFF	220>420 mm circumference	
HEIGHT	80 mm	
WIDTH	103 mm	
DEPTH	129 mm	
WEIGHT	250 g without batteries	
ACCURACY	Blood Pressure ±3 mmHg	
	Pulse Rate ±5%	
CONFORMS TO	ARTG 137264	
WARRANTY	1 year	